


MPK Protectors and Cosmesis

Protection for everyday life



Functional Cosmesis for C-Leg 4[®] and Genium[®]

The Functional Cosmesis consists of a durable functional knee part, a functional foam shank that can be individually shaped, and an attractive, easy to care for functional stocking. The combination of a great look and functionality will not only please your patients, it also makes your fabrication work easier.

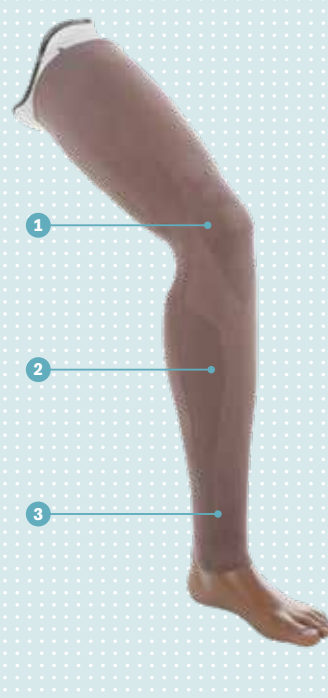

Functional Knee Part

- Exact fit and robust
- Significantly reduced effect on knee function compared with conventional standard solution
- Spares the difficult process of fabricating a functional knee solution in-house
- Forms a natural-looking knee part
- No conventional repair work such as gluing tears with one-piece foam solutions

Functional Shank

- The inner contour is designed precisely to fit the shape of the knee joint
- Makes donning and doffing considerably easier
- Firm fit at the knee joint makes sanding easier
- Natural feel

Functional Stocking

- Washable
- Functional zones: elastic knee part, concealing fleece yarn zones
- Natural shading
- Splash and dirt repellent

Your benefits at a glance:

- High degree of prefabrication
- Highly durable
- Improved hygiene and protection
- Natural look

Advantages for your patients:

- Natural look and function combined in one solution
- Inconspicuous appearance in long or short clothing
- Highly functional in everyday life

Order Information

Select Article Numbers

Product	Article Numbers
C-Leg 4	3F1=1 + 99B120=*
Genium	3F1=2 + 99B120=*

Select Stocking Size


MTP to floor measurement	Socket circumference (below ischium)	
	15.7-18.1 in	18.1-22 in
14.9-18.9 in	S	L ¹
>18.9 in	L	L

¹These combinations of measurements can occur in prosthesis users who are very tall and also very slim. For the Functional Stocking to fit the socket well, the socket diameter must be enlarged by adding material in these cases.







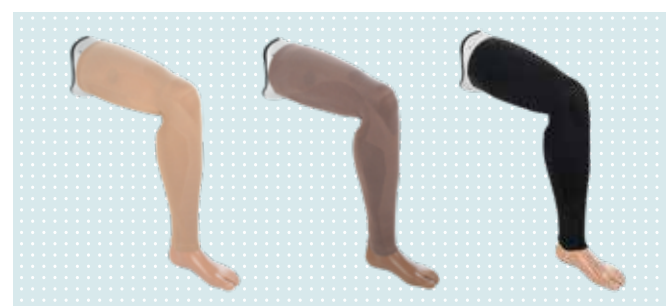

Select Color Code

beige = 4
brown = 15
black = 7

Order example for C-Leg 4: 3F1=1 + 99B120=L-4



Cosmetic Foam Cover and Protector Cover Options

	Kenevo	C-Leg 4	Genium	X3
Cosmetic Foam Cover	 3S26=R44 3S26=L44	 3S26=R44 3S26=L44 3F1=1	 3S26=R44 3S26=L44 3F1=2	
Protective Cover	 4X840	 4X860=S 4X860=M 4X860=L 4X860=S-8.4 4X860=M-8.4 4X860=L-8.4 4P862	 4X880=S 4X880=M 4X880=L	 4X193-1 4X900
Shield for Protective Cover		 4P863=1 4P863=2 4P863=3		

