

Active TH Harness 21A47

The 21A47 soft harness is intended for transhumeral fittings with myoelectric or passive elbows. The support, which consists of an arm sling and a sleeve, is completely removable and is not fastened to the socket with a strap. Users can put it on and take it off with one hand. Soft padding in the armpit and breathable materials, which are fully washable, help make the support comfortable.

The belt guide at the back allows the arms to swing naturally, supporting more physiological movement. A hook-and-loop fastener in the front enables the user to make fine adjustments in certain situations, for example, when carrying something heavy. The arm loop (axilla cuff) portion of our Active TH (trans-humeral) soft harness is also available separately. This provides users with the flexibility to mix and match sizes as well as the option to only replace the arm-loop as necessary.



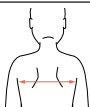
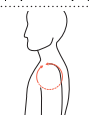
Features and user benefits of the Active TH Harness:

- Textile materials
 - High level of comfort
 - Moisture-wicking
 - Machine washable
- Axillary padding in the arm pit
 - Enhanced comfort
 - Improved pressure distribution
- Cord guide on back
 - Allows physiological movement of contralateral arm
 - No fixation to the socket
 - User can put on and take off the harness with one hand
- Adjustment buckle on front side
 - Users can tighten or loosen harness when necessary. For example, when carrying heavy objects.
 - May be modified for additional control straps, etc.
- Versatility in installation/usage
 - Users can mix and match sizes
 - Arm Loop is available separately. Allowing replacement of only the arm-loop as necessary

Technical data

| Reference number | 21A47 | 21A48 |
|------------------|--|---------------|
| Description | Transhumeral soft harness for electronic or passive elbows | Arm Loop Only |
| Size | S, M, L | S, M, L |
| Versions | Left, Right | Left, Right |
| Color | Grey (1) | Grey (1) |

Size chart

| | Side | | Size | | |
|---|------|------------------------------|-----------|-------------|-----------|
| | | | S | M | L |
|  | L/R | Back width | 30–36.6cm | 36.6–43.2cm | 43.2–50cm |
|  | L/R | Armhole circumference | 38–45.3cm | 45.3–52.6cm | 52.6–60cm |

Order Examples

| Article Number | Description | Side | Size | Color |
|----------------|-------------------|-------|--------|-------|
| 21A47=L-S-1 | Active TH Harness | Left | Small | Grey |
| 21A48=R-M-1 | Arm loop only | Right | Medium | Grey |